neno[®]

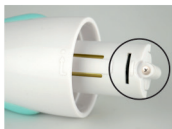
DENTI

Table of contents

Spis treści

- 4** Instrukcja obsługi
- 7** User manual
- 10** Bedienungsanleitung
- 13** Návod k použití
- 16** Használati utasítás
- 20** Používateľská príručka
- 23** Vartotojo vadovas
- 26** Lietotāja rokasgrāmata
- 29** Kasutusjuhend
- 32** Manual de utilizare
- 35** Brukermanual
- 38** Brugermanual
- 41** Käyttäjän käsikirja
- 44** Användarhandbok

A**B****D****C**

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EN

USER MANUAL

Dear Customer, thank you for purchasing NENO DENTI. The product you have obtained is an electric toothbrush for small children. This device is used to brush the teeth of a baby and maintain their oral hygiene. Before using the device, please read instructions below.

01. PRECAUTIONS

1. If you do not understand the instructions for use or precautions listed in this manual, do not use the device. Incorrect use of the device may damage it.
2. Do not disassemble the product. This may result in damage to the product, electric shock and voiding the warranty.
3. Do not use accessories other than those recommended by the manufacturer.
4. Do not use flammable or explosive substances to clean the product.
5. This product is only used for brushing teeth and massaging gums. Do not use the device for any other purpose.
6. Stop using the device if it is clearly damaged.
7. Do not submerge the product in water. This can lead to a malfunction and may result in fire, electric shock or explosion.
8. Keep the device away from water containers such as a bath or sink.
9. Do not replace the battery of the device when your hands are wet.
10. Do not use used, old batteries in your device.
11. Pay attention to the poles of the battery when placing it in the device.
12. Do not place the device in rooms with high moisture or high temperatures, such as bathrooms.
13. Remove the battery from the device if you will not use it for an extended period of time to avoid the risk of fire.

14. If the toothpaste used contains peroxides or bicarbonate (often present in tooth whitening pastes), wash the device thoroughly after use to avoid destroying its plastic casing.
15. If using the device causes your child discomfort or bleeding from the gums, consult your dentist before continuing to use it.
16. Children may only use the appliance under supervision adults who have read the operating instructions.
17. The device is not a toy – keep it away from children, so that they do not use it without the supervision of an adult.
18. The device is used for personal hygiene and should not be used on many people in dental offices.

02. DESCRIPTION OF THE DEVICE

SEE FIGURE A

1. Head with bristles
2. Light
3. Power switch
4. Screw cap/ battery space
5. Suction cup

03. PACKAGE CONTENTS

SEE FIGURE B

1. Toothbrush
2. Head with bristles
3. User Manual
4. Massaging heads:
 - 4.1. Rubber head for gum massage
 - 4.2. Rubber head for tongue
 - 4.3. Rubber head for cleaning gaps between teeth

04. POWER

The device needs 1 piece of AAA battery for proper operation. Unscrew the lid at the bottom of the device to access the battery space. The battery compartment cover is held in place with a screw so that the child does not have access to the power supply of the device. Unscrew the cover with a screwdriver and put in the AAA battery as shown in the diagram. **SEE FIGURE C**

Place the lid back in place and turn it, making sure it stays in place. The bump on the lid should be under the bottom of the arrow on the housing before turning, and under the dot after turning the lid. **SEE FIGURE D**

05. USING THE DEVICE

1. Place the brush head onto the device. It should be put on the plastic part with the light inside it.
2. Press the power button of the brush to turn on the light of the device head.
3. Squeeze toothpaste onto the brush.
4. Press the toothbrush gently to your teeth and gums.
5. Press the power button again to turn on the brush vibration and start thoroughly cleaning the baby's teeth. Every 30 seconds of brushing, the toothbrush stops vibrating for 1 second to allow the user to Change the brushed space. Correct brushing of teeth should take a total of 2 minutes. The device will turn off automatically after 2 minutes.
6. When you're done brushing your teeth, turn off the toothbrush by pressing the power button again.
7. After cleaning the baby's teeth, wash the toothbrush under running water. Remove the remaining water from the bristles of the brush and use soft cloth to wipe the handle of the device.
8. When you are not using the brush, keep it upright.

NOTE: to avoid toothpaste splashing around, turn off the device before removing it from the baby's mouth For the sake of your child's hygiene and safety, replace the brush head every three months.

06. CLEANING THE BRUSH

1. Wash the brush head and its bristles after using the device and let it dry naturally.
2. Do not press the bristles or pull them out with your fingers. This can deform the toothbrush and destroy it.
3. Wipe the handle with a soft piece of cloth and let it dry naturally.
4. Do not submerge the handle underwater.
5. Do not use water above 40 °C or acidic cleaning agents to clean the device.

07. FAQ

Q: Vibrations weaken during use:

A: You push the toothbrush too hard. It should gently touch the teeth.

08. SPECIFICATION

Power source: 1x AAA battery

Waterproof: IPX7

Nylon bristles: Yes

Strokes per minute: 17000

Built-in timer: Yes

Built-in LED light: Yes

Dimensions: 149x28x28

Weight: 28g

09. WARRANTY CARD

The product is covered by a 24-month warranty. The terms of the guarantee can be found at: <https://nenopl/warranty>

Details, contact and website address can be found at: <https://nenopl/contact>
Specifications and contents of the kit are subject to change without notice. We apologise for any inconvenience.

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EN

The crossed out trash can symbol indicates that unusable electrical or electronic devices, its accessories (such as power supplies, cords) or components (for example batteries, if included) cannot be disposed of alongside with household waste. In order to dispose of the devices or its components (for example, batteries) deliver the device to the collection point, where it will be accepted free of charge. Disposal is subject to the recast version of the WEEE Directive (2012/19/ EU) and the Directive on batteries and accumulators (2006/66 / EC). Proper disposal of the device prevents degradation of the natural environment. Information about the collection points of the facilities is issued by the competent local authorities. Incorrect disposal of waste is subject to penalties provided for by the law in force in the given area.

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